


VISIONA

connected ultra-low hotel-style medical bed,
REF. COV58439BTD

AHF

A MEDJUNIC COMPANY

The **VISIONA** ultra-low electric medical bed is a "connected" bed featuring 5 types of alarm
> 1 main alarm (included in the package): **bed height - indicates incorrect bed height adjustment to limit the risk of resident falls and injuries**
> 4 optional alarms: resident presence, bed rail position, brakes off and incontinence detection.
Intuitive operation thanks to the centralised controls on the Nurse Control unit fixed onto the bed.
The sensors necessary to resident safety are invisible, so there is no disruption in the room environment.

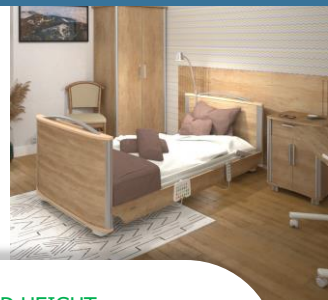

Technical specifications

HEIGHT ADJUSTMENT

LIFTING COLUMN < SPECIFIC TO FLOORE, VISIONA >

• The design with electric actuators integrated into the columns allows totally free access underneath the bed, improving the work of the staff and enabling productivity gains.

FUNCTIONS

MULTIFUNCTIONAL

• Variable height / Backrest (with autoregression) / Footrest (knee break) / TR/ATR / Cardiac chair position



ALERT – BED HEIGHT

(INCLUDED AS STANDARD)

• Detects if the bed is not in the safe low position (23 cm from the floor)

Green light: bed in safe low position

Red light: bed in an unsafe position

The visual alarm can be accompanied by an audible alarm

HEAD AND FOOTBOARDS

- Customisable: wide selection of laminates and varnished solid woods

ACCESSORY ANCHORS

- 4 for drip holder
- 2 for hand control pendant holder (at the head of the bed)
- 4 for trapeze bar holder (at the head & foot of the bed)

WALL BUMPERS

- Aluminium wall bumpers (supplied) are easily screwed into the headboard

CPR

- The emergency bed flat function can be quickly activated using 2 easily accessible bilateral handles

BED BASE

- 5 removable slats made of recycled ABS
- (optional) Compact laminate flat bed base

(optional) SIDE RAILS

- 3 concepts: half-length rails, 3/4-length rails, or full-length rails.

CONTROLS

- Intuitive wired hand control pendant (included), with the option to disable certain functions
- Nurse Control unit for caregivers (included)
- The edge flashes red if an alarm goes off

EXTENSION

- Integrated in the foot of the bed with mattress extension rest (easy access and tool-free setup)

BRAKES

- Semi-centralised, bilateral at the head/foot of the bed

CASTORS

- Swivel castors (360°), dual wheel, diam. 50 mm
- Encased as standard for a "non-medical" look

FRAME

Space under the bed totally clear < SPECIFIC TO FLOORE, VISIONA >

- The height adjustment system allows completely unobstructed access underneath the bed

Key data

Variable height	23 cm – 68 cm
SWL (patient weight)	220 kg (185 kg)
Empty weight (without accessories)	100 kg
Backrest angle	70° (max.)
Knee break	30° (max.)
TR/ANTI-TR angle	17° / 17° (max.)
Extension	12 cm
Dimensions (without side rails)	220 x 106 cm
Lying surface	200 x 90 cm
Materials	Steel, Solid wood, Laminate
Dual wheel castors	360° / Ø 50 mm
Supply voltage	100-230 V, 50/60 Hz Max. 3.15 A
Certification Standards	CE (REG (EU) 2017/745) NF EN 60601-2-52 NF EN 60601-1-2 NF EN 60601-1

Customisation and Accessories



Vicenza Light Oak
<NEW> MARCH 2024



Nebraska Oak
Natural



Nebraska Oak
Grey

See our colour chart for all the laminate colours available

Colour (bed frame)



Cream white
RAL 9002



Anthracic grey
RAL 7015
FROM 2025

Side rails available

- Half-rails
- 3/4-length rails
- Full-length rails

Hygiene & Disinfection

- Cleaning products must be sprayed on or applied with a disposable cloth
- It is also recommended to start cleaning at the top and work down

Accessories available

- Trapeze bar with retractable strap
- 4-hook drip pole
- Flexible pendant holder
- Under-bed night light
- Wall bumper, length 19 cm
- Battery
- Urine bag holder (on 3/4-rail)
- Grip handle
- Compact laminate flat bed base
- Mattress pump bracket
- Standard and large loop strap attachments
- Reading lamp
- Sound kit
- Pillow net
- Foam side rail pad
- Extension for 3/4-rail
- Side panel trim for 3/4-rail
- Wood trim for 3/4-rail
- Transport trolley

Optional – 4 additional alarms



ALARM – SIDE RAIL POSITION

OPTION 1

- Side rails are essential to ensure the safety of residents, especially during the night. The alarm goes off if any of the side rails are in the low position.

OPTION AVAILABLE FOR HALF-RAILS & 3/4-RAILS ONLY

REF: LSUGHCPB3/4LF – Side rail detection

REF: LSUGH2BPLTXO1/2G2C

(Set of 4 connected half-rails)

REF: LSUGH2BPLTXO6C

(Set of 4 connected 3/4-rails)



ALARM – BRAKES OFF

OPTION 2

- The bed brakes being off can be an aggravating factor in a fall or trauma injury. There are sensors on each axle and the alarm goes off if the brakes are not on.

REF: LSUGHCLF – Brake detection



ALARM – RESIDENT PRESENCE

OPTION 3

- If a person is confused or disoriented, it is important to be able to detect whether they are in bed or not. Depending on the resident profile, this alarm can be customised to go off after 15 seconds, 3 minutes and 12 minutes.

REF: LSUGHCPPF – Presence detection



ALARM – INCONTINENCE DETECTION

OPTION 4

- Incontinence affects many disoriented persons and can be difficult to detect. Our system alerts care staff as soon as traces of moisture are detected in the bed.

REF: ALARME – Incontinence detection

REF: 40116600 – Incontinence mat



3 VISUAL ALARM INDICATORS

- VI Box flashes red
- Nurse Control flashes red round the edges
- LED for the function on the Nurse Control flashes

Focus – Nurse Control

Bed functions

- The traditional bed functions are available on the panel and be individually locked

Additional positions

- CPR emergency position
- Trendelenburg emergency position
- Cardiac chair position
- Night position (low position)



You can choose to use the audible alarm or not

- To enable or disable it, press the corresponding button. The audible alarm is enabled when the LED is on

Detection alarm

- Each alarm is represented by a green button. To enable or disable an alarm, press the corresponding button. The alarm is enabled when the LED is on



Want to know more?

- Watch our specially made video to understand and visualise the different alarms offered on the VISIONA connected bed.
Click on the icon to find out more